


Easy Plate™

Save time, labor, and cost with Kikkoman Easy Plate.

Easy Plate is a convenient, prepared media plate for faster and easier microbiological testing.

Easy-to-use,
does not require
a plate spreader.

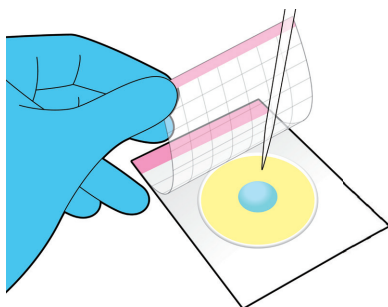

Product	Easy Plate AC Aerobic Bacteria
Incubation	35±1°C, 48±2 hours
Shelf life	18 months after manufacturing date
Certification	AOAC PTM, MicroVal, NordVal



Product	Easy Plate CC Coliform
Incubation	35±1°C, 24±1 hours
Shelf life	18 months after manufacturing date
Certification	AOAC PTM, MicroVal, NordVal



Product	Easy Plate EC E. coli & Coliform
Incubation	35±1°C, 24±1 hours
Shelf life	18 months after manufacturing date
Certification	AOAC PTM, MicroVal, NordVal



Product	Easy Plate SA Staphylococcus aureus
Incubation	35 or 37±1°C, 24±1 hours
Shelf life	12 months after manufacturing date
Certification	AOAC PTM, MicroVal, NordVal



Product	Easy Plate YM-R Yeast and Mold rapid
Incubation	25°C±1, 48±2 hours
Shelf life	18 months after manufacturing date
Certification	AOAC PTM, MicroVal, NordVal



Product	Easy Plate EB Enterobacteriaceae
Incubation	35±1°C or 37±1°C, 24±1 hours
Shelf life	18 months after manufacturing date



Product	Easy Plate AC-R Aerobic Bacteria Rapid
Incubation	35±1°C, 24±1 hours
Shelf life	12 months after manufacturing date

Colony Counting System for Easy Plate

Download for free

Automatically counts the colonies in imaged Easy Plate. Simpler and time-saving way to count colonies than manual counting advances efficiency and eliminates the effects of variation and skill level among operators.

Download software from here.



Boost lab productivity easily with automatic colony counter

Learn more about software from here.



The screenshot displays the CFU software interface. On the left, a data table shows 'Counted date: 2023/09/05 17:32:01', 'S/N: AC211207-01-010828', 'Filename: easyplate_ac_02', 'Type: AC', and 'Status: OK'. The center features a petri dish image with red square markers indicating colonies. On the right, an 'Edit Mode' panel shows 'CFU' selected, 'AC' as the type, and '158' as the count, with 'Save' and 'Clear' buttons. Below this is a 'Show marks' toggle set to 'AC'.

- Up to 7 sheets can be inserted at once.
- Highly accurate analysis with AI technologies.
- Automatically determines which of the Easy Plate products is used and measures colonies.
- It also identifies spread colonies.
- Parameters to judge colonies are available with no initial settings.
- Able to modify the result and add the comments.

Supported Languages:

English, German, Spanish, French, Korean, Thai, Vietnamese, Simplified Chinese, Traditional Chinese, Japanese



*It requires the scanner (ADS-4300N by BROTHER INDUSTRIES, LTD.) for imaging Easy Plate. We do not sale the scanner.

kikkoman

Kikkoman Biochemifa Company

2-1-1 Nishi-shinbashi, Minato-ku,
Tokyo 105-0003 JAPAN
biochemifa@mail.kikkoman.co.jp
<https://biochemifa.kikkoman.com/e/>

© 2025 Kikkoman Corp. (6012Y250201)

Follow us on
social media!

